



MIICRAFT

NEW!



THE LEADING
3D PRINTER
FOR DENTISTRY

COST EFFECTIVE
TRUE OPEN SOURCE
INDUSTRY LEADING SUPPORT

**Miicraft 125
Ultra DLP printer**

See Page 7

bredent group



red & white with crea.lign®

See Page 9

keystone
industries



vhf

NEW!

vhf milling and grinding mills
"Dental from Deutschland"

See Page 6



NEW!

**Enamelite Low Fusing
Ceramic Spray Glaze**



See Page 12
(Back Page)

QUALITY • RELIABLE • ECONOMICAL



MESTRA Model Trimmer

It features a powerful 1600 W engine with integrated brake, opening lid for easy access to the disc, operation warning light, disk and more...

SPECIAL
\$895⁰⁰
080093



easily
removable
cover



• velcro disc
• easy to replace

IRIS 2 Vacuum Mixers

IRIS 2 Model 080516

Simple operation mixer just place bowl & select mix time.

PRICE
\$1195⁰⁰
080516



IRIS
Evolution

Programmable

PRICE
\$1295
080517



IRIS
Advanced

Programmable

PRICE
\$1395
080518



Polymerizer M-9

- Holds up to 9 flasks
- 2 stage programmable



PRICE
\$2495
080400

Polymerizer M-18

- Holds up to 18 flasks
- 2 stage programmable



PRICE
\$3195
080402

Polymerizer E2

- Ready to shower at all times
- Holds up to 18 flasks



PRICE
\$5440
080600

MORE PRODUCT SPECIALS

BDM Knives

#6A & #7A autoclavable knives feature tough plastic handles and stain resistant blades for sterilization durability.



#65570 (6A) **\$24.25**

#65600 (7A) **\$22.65**

AUTOCLAVABLE INSTRUMENTS

BDM Spatula

Autoclavable Spatula



\$24.25

#79480

C/D Fahnstock Wax Knife

Autoclavable all metal wax knife



\$46.95

#HSL-109-17

PLASDENT Mixing Bowls

Spectrum flowbowl mixing bowls 600 cc large.

\$6.95
Red

#904MBL-2C Blue transparent

#904MBL-5C Red transparent



\$6.95

Blue

**BUY 3 BOWLS RECEIVE
FREE NYLON* SPATULA!**

*Nylon Spatula



\$3.95

#904SA

BDM Plaster Traps

Trap-Eze Disposable Plaster Trap

- New smaller size (9.5") height for smaller spaces.
- Sealed bucket elements odours
- Easy to change and install

NEW SMALLER SIZE!

#62170 (2 gallon complete)

\$59.95

#62180 (2 gallon refill)

\$26.95



C/D Labo-Mate 80

Articulator

- Magnetic split cast insures easy mounting and removal of models.



\$34.95

#02030

SAVE \$10

Reg. \$44.95

C/D Labo-Mate 150

Plasterless Articulator



\$79.95

#02060-SET

SAVE \$30

Reg. \$109.95

TOKUYAMA MR. BOND

- Primer for bonding self-curing acrylic resins to non-precious metals.
- Strong adhesion

\$90.45

#20133



TOKUYAMA Rebase II

Chairside Hard denture reline material

- Less heat generation and irritation than competing products.



\$23^{value}
#5410.060

**Get a FREE Bur
with Purchase!**

#20506

\$180.90

Sunburst Micro-Motor

- 30,000 rpm light duty micro-motor
- Compact size



SPECIAL \$350⁰⁰
Reg. \$425.00
#P-900/X-99L

BUY 2 Receive an Extra 12% OFF

DIGITAL WORKFLOW SOLUTIONS FOR DIGITAL DENTISTRY

EXTRA-ORAL SCANNING



open technologies

Optical 3D Scanner

Escape the tyranny of annual license and maintenance fees!

Designed and made in Italy, Open Technologies offers dental professionals a range of highly precise and robust extra-oral scanners to choose from and options that allow you to spend only what you need to.

All Open Technologies scanners come with ZERO annual hardware or software license, maintenance fees and offer seamless integration with exocad DentalCAD.

*all pricing excludes PC and CAD software.

eaSy

Perfect as an add-on scanner to permit for extra capacity or as a cost-effective starter unit where the plan is to focus strictly on crown, coping and bridge manufacturing.



Starting from

\$12,400

smaRt Small or BIG

Add precision and extra standard features for those with more complex scanning requirements such as implant work. Opt for the Smart Big where being able to scan a full-size articulator is a priority.



Starting from

\$14,500

deluXe

The best value offering. All the features of the Smart Big plus colour scanning and extra standard modules.



Starting from

\$16,500

neWay

Fastest scanning speeds and unique design, coupled with the most recent software interface.



Starting from

\$18,750

DIGITAL WORKFLOW SOLUTIONS FOR DIGITAL DENTISTRY

CAD DESIGNING

exocad

Dental CAD

High tech tools for high precision artists!

EASY TO LEARN


The strength of exocad DentalCAD is its simplicity. It's extremely logical and easy to use, even for first-time users. The software emulates the analog work of a dental technician.

This will help you minimize training costs and maximize productivity and profits. Exocad DentalCAD is reliable and robust even when dealing with complex cases on a daily basis.



EASY TO UPGRADE

Expand your exocad DentalCAD standard version according to your needs. Exocad add-on modules enhance the feature-set of the standard version. This modular system allows great flexibility and you can purchase (and spend) just what you need.

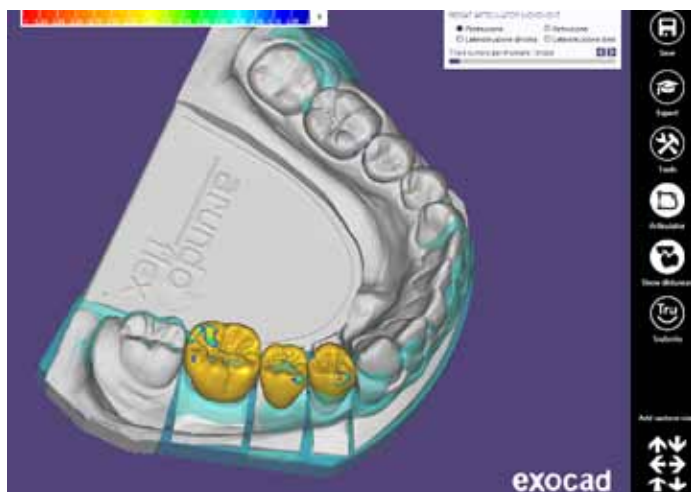

DIGITAL WORKFLOW SOLUTIONS FOR DIGITAL DENTISTRY

MILLING AND GRINDING

vhf

VHF milling and grinding mills - "Dental from Deutschland"

VHF Camfacture AG has been making milling and grinding machines under private-label for companies around the globe and under it's own brand for thirty years. VHF dental milling and grinding machines are robust workhorses, designed to reliably withstand the stresses created during milling and grinding; especially of harder materials like chrome-cobalt, glass ceramics and PMMA.

VHF offers a fully integrated solution with validated work-flows. Companies in the VHF group manufacture all components of the machines including the internal electronics, the milling and grinding tools, and the dental CAM that define the milling and grinding strategies. This ensures that outputs from the machines are exactly as expected and desired, making manufacturing processes predictable.

K5 - Dry milling of zirconia, chrome-cobalt, PMMA, wax and polymers

The K5 is a 5-axis machine designed for precision and longevity. A massive bed of cast aluminum means the K5 is extremely stable during milling and as such is one of the few bench top machines robust enough to handle milling of chrome-cobalt. A special design creates an internal vacuum in the milling chamber, drawing damaging debris and powders away from the spindle, contributing to extended spindle life. The K5 b-axis will rotate up to 35 degrees, making the milling of thicker elements like models a snap. The higher rotation range also allows for the milling of discs up to 40mm thick.



Starting from
including DentalCAM7 **\$43,500**



N4 - Wet Grinding of glass ceramics and titanium

The N4 is a compact 4-axis wet grinding machine suited for the production of glass ceramics, zirconia and titanium. The N4 has a fully integrated liquid cooling system, simplifying the overall equipment needs and setup.

\$35,735 Starting from
including DentalCAM7

Dental Direkt zirconia - German zirconia for the German mill

BioZX2 from Dental Direkt is the perfect accompaniment for VHF mills, or for any mill as a matter of fact. BioZX2 is a high-translucent, high strength material that is isostatically pressed, resulting in a more dense material that mills more predictably without chipping. BioZX2 comes in white discs that are dipped/shaded, or in pre-shaded discs based on the Vita Classical shade system - various disc diameters are available.



\$110 Starting from
10mm discs

DD BioZX²

**Dental
Direkt**

DIGITAL WORKFLOW SOLUTIONS FOR DIGITAL DENTISTRY

RAPID PROTOTYPING

MIICRAFT

Miicraft 125 Ultra DLP printer

The latest DLP printer from Miicraft is the perfect upgrade unit from inexpensive "starter printers", or simply as the initial entry point for someone making that first move into rapid prototyping for small scale custom dental production. The heart of the Miicraft is the powerful and superior quality light engine, which operates at far below maximum output, resulting in prolonged service life.

Capable of printing virtually any dental CAD design, the Miicraft 125 is a fully open and adjustable printer that will work with literally any print resin on the market that cures at 405nm. As new printable materials come to market your Miicraft printer will always be open to them, effectively future-proofing your investment.



**THE LEADING
3D PRINTER
FOR DENTISTRY**

BRING DIGITAL DENTISTRY
TO YOUR CLINIC TO INCREASE
PRODUCTIVITY & REDUCE COSTS

COST EFFECTIVE
TRUE OPEN SOURCE
INDUSTRY LEADING SUPPORT

Starting from
US Dollars **\$8,995**

CADWorks nesting software – designed specially for the Miicraft 125 Ultra

CADWorks nesting software simplifies and speeds up the preparation of design files. Easy and intuitive, CADWorks software allows users to prepare supports using powerful automatic tools, or to do so manually when doing so is more effective.

CADWorks nesting software is included with the Miicraft 125 or available for sale starting at \$1,000 per licence.



POWERED BY
CADworks 3D

Keyprint print resins – “Juice for Your Printer”

Now available for sale in **Canada**

Need a print resin for your open printer? Keyprint resins from Keystone Industries will do the trick.




Keymodel beige - **\$325** per kg.

Keyortho model grey - **\$325** per kg.

Keymask pink : soft tissue model resin - **\$285** per kg.

Keycast: casting resin - **\$535** per kg.

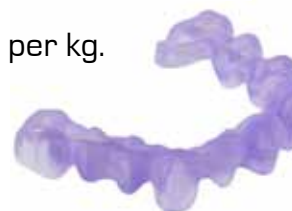

keystone
industries

keyprint

Coming soon: Keyguide clear surgical guide material, Keysplint hard splint material and Keysplint soft splint material.



neo.lign teeth

- Natural esthetics, following nature's example
- Anterior and posterior teeth are constructed in thin layers from high impact PMMA Composite (polymethyl methacrylate).
- Available in classic A-D range and BL3 bleach
- Neo.lign posterior teeth enable implementation of all occlusion concepts due to the multi-functional chewing surface of G and L design.
- With the neo.lign L designs, the viso.lign system provides a "lingualised occlusion design" specifically for manufacturing full dentures.

Special working mould guide offer only \$595.00

Zi-finish

- Pre-polishers for smoothing and pre-polishing; available in 3 shapes- lens, wheel and pointed cone
- Polishers for high lustre polishing also available in 3 shapes- lens, wheel and pointed cone
- Diagen turbo green stones

save \$33

Lens	Ø 20 mm x 1.25 mm
Wheel	Ø 11 mm x 2.3 mm
Cone, pointed	Ø 4.2 mm x 12 mm

Special
\$185
Reg. \$218

Advantages of the new Zi-finish kit:

- Enables fast surface polishing of zirconia
- 2 stage system of pre-polish and polish- simplifies the work and reduces work time.
- Lasts longer than others on the market
- Can also be used with ceramic and non-precious metals



Kit
#33000836



red & white with crea.lign®

Red & white with crea.lign® is the solution for:

- Implant-supported, fixed or removable permanent restorations
- Natural customisation of full dentures
- All substructure materials such as BioHPP® and PEEK
- Complete coating even of full-arch restorations (e.g. Prettau bridges)
- Can also be used with ceramic and non-precious metals

SPECIAL
\$595⁰⁰

Reg. \$675.00

CLIGNSETG

Diatit Tungsten Carbide Burs

Diatit tools feature the advanced Diatit wear protection coating. This is a special material which is added into the surface of the bur after it has been produced. It hardens the tool surface and reduces the surface friction.

Ask Us About Our Catalog For A Full Range Of Burs Available!

**Buy 5
Get 1
FREE!**

Of equal or
lesser value*



\$28⁹⁵

REF D200 KG 23



\$48⁵⁰

REF D237 KG 65



\$50⁹⁵

REF D263 KG 60



\$53²⁵

REF D274 KG 60



\$28⁹⁵

REF D292 KG 23

**Smooth
Running
& Longer
Lasting**

INFECTION CONTROL SPECIALS

FLIGHT Clave16

- Fully automated 16 L chamber with pre-set programs.
- Stainless steel chamber
- External printer is an option
- Internal memory of last 20 cycles



SPECIAL
\$4695⁰⁰
Reg.\$4950

A9001-16

FLIGHT Clave23

- Fully automatic 23 L chamber with pre-set programs.
- Stainless steel chamber
- External printer is an option
- Internal memory of last 20 cycles



SPECIAL
\$5650⁰⁰
Reg.\$5950

A9002-23

FLIGHT ClaveB

- Class B autoclave with pre/post vacuum
- Fully automatic 23 L chamber with pre-set programs.
- External printer is an option



SPECIAL
\$8360⁰⁰
Reg.\$8800

A9003-23B

Tuttnauer EZ9/EZ10

Autoclaves will meet all your sterilization needs

- Chamber door closed during entire cycle
- Hepa filtered air is circulated in chamber during drying
- All program parameters can be customized and stored
- 2 year parts and labour warranty
- Easy to read panel



EZ9 - PRICE
\$7330⁰⁰

EZ10 - PRICE
\$8295⁰⁰

Tuttnauer Valuklave 1730

- Compact size[7" chamber]
- Fast cycle time - (11 minutes from a hot start)
- Automatic shut-off at the end of both the sterilization and drying cycles.
- 1 year parts and labour warranty



#1730
PRICE
\$3750⁰⁰

Nitrile Powder-Free Exam Gloves

CODE/SIZE
92395-XS
92396-S
92397-M
92398-L
92399-XL



200
Gloves / Box

10 + 1

Purchase 10, Get 1 FREE!

*Double the
Normal Amount! (100/box).*

Earloop Face Masks

5 + 1

**Purchase 5,
Get 1 FREE!**



Blue
(50/box)

#5614

CROSSTEX Confirm 10

- In office biological monitoring system provides final results in just 10 hours
- Simply process test indicator vial in a normal steam sterilization cycle, remove, crush ampoule and place in incubator along with an unprocessed control vial for required 10 hour incubation period.




CROSSTEX Sure-Check® Strips

- For use in steam and chemvapor
- Multi-parameter indicates if all conditions for sterilization have been met, time, temperature and presence of steam or chemical vapor.
- Can be used in pouches, tubing, autoclave bags



CROSSTEX Sure-Check® Self-Seal Sterilization Pouches

Built-in multi-parameter internal/external process indicators that change color when all three criteria for sterilization have been met – time, temperature and pressure [steam].

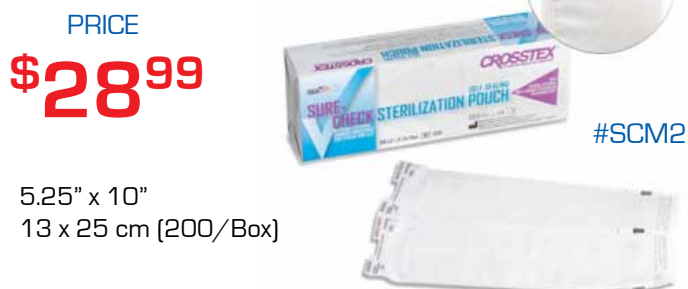

CROSSTEX Sure-Check® Sterilization Pouches (Steam)

Sure-Check® Sterilization Pouches feature internal/external multi-variable indicator technology that confirms all three criteria for sterilization have been met.



MICRYLIUM BioTEXT®

BioTEXT is the only clinical disinfectant specifically designed for use on delicate surfaces. It is proven to kill benchmark organisms on naugahyde, leather, vinyl and rubber surfaces. Ideal for tubing, fabrics, carpets, touch screens and soft plastic surfaces such as keyboards.

LIQUID **\$69⁹⁵**

- One step cleaning power
- Gentle formula for use on delicate and textured surfaces
- Complete stain and odor elimination

01-TXT-005 (5L Box)

WIPES *NEW* **\$116⁹⁹**

- Biodegradable and absorbent
 - Strong non-woven material designed for use on irregular surfaces
 - Resealable package keeps product fresh
- 04-TXWI-802 (8 Rolls & 2 Canisters)



MICRYLIUM BioSURF®

BioSURF is the world's fastest one step surface disinfectant. Unlike most disinfectants that can take up to 10 minutes to disinfect, BioSURF works within seconds. Despite the impressive performance, BioSURF is designed to not harm people or the environment.

PRICE **\$92⁹⁵**

01-SURF-005
(5L Box)



ENAMELITE Low-Fusing Ceramic Spray Glaze

Glaze ceramics like never before!

- Reduces labor cost, increases efficiency
- Spray a thin, uniform, & accurate coat IN SECONDS
- Faster drying than brush-on applications
- Compatible with ALL dental ceramics

keystone
industries

PRICE

\$115⁰⁰

#5100179



Aquasil Ultra

The perfect restoration always begins with an Aquasil Ultra Impression!



Dentsply
Sirona

4+1

Purchase 4, Get 1 FREE!

NEW!

AeraMax PRO AM III - Stainless Commercial Air Purifier

- True HEPA filtration captures 99.97% of airborne contaminants including viruses and allergens
- Activated carbon filter reduces odors
- Patented EnviroSmart technology detects the environment and automatically adjusts performance.
- Cleans the air of 300 to 700 square feet spaces
- 5 Year Limited Warranty
- Floor or Wall mounted

#9450001 (with stand)

\$1169⁰⁰

#9416201 (wallmount)

\$1119⁰⁰

