

Spring 2026

# Flyer

Promotions in effect April 1st - July 31st 2026



Scan to place your order or speak to a representative or call toll-free:

1-800-268-4442

**C/D** *Central Dental*  
The laboratory and denturist supply specialist.

# Endless possibilities

Blue Line®



## SR Vivodent® DCL

- 42 individual anterior tooth moulds
- Made from durable, highly cross-linked DCL material
- Contemporary shade range



LEARN  
MORE

# Precision at the push of a button

## IvoBase<sup>®</sup> System

Denture base material with fully automated injection process

- High precision due to shrinkage compensation
- Fully automated process
- Predosed units for an accurate mixing ratio

©2026 Ivoclar Vivadent, Inc. Ivoclar and IvoBase are trademarks of Ivoclar Vivadent, Inc.  
\* Promo Code: 722 IvoBase product code 628518BU

ivoclar.com  
Making People Smile

Save 20%  
on an IvoBase<sup>®</sup>  
(628518BU)\*

ivoclar

**Renfert**

making work easy

# Renfert Essentials: Core Solutions for Modern Dental Labs



## MT3 Marathon

(Partially diamond-coated disc)

Item No.: 18081500

~~\$3,204.00~~

**\$2,895.00**

SAVE OVER \$300



## MT3 Klettfix

(80 grit, Corundum disc)

Item No.: 18081000

**\$1,950.00**

**+ FREE**

## Klettfix Disc 5 pk

With purchase of MT3 Klettfix

Item No.: 67000100



## SYMPRO

Automatic  
Appliance Cleaner

Item No.: 67001000

**\$1,195.00**

**+ FREE**

**20 Cleanings!**

help:ex plaque f

With purchase of SYMPRO

Item No.: 67000100 (4 x 1 Liter)



## Renfert Firing Paste

Firing paste for stabilizing ceramic works

Item No.: 19220304 (4 x 12g)

Renfert's new Firing Paste is the smart choice for dental technicians who demand consistency, clean results, and zero compromise — even with the most delicate ceramics. Usable up to 1,500°C.

**\$52.95** (4 pk)



Dynex  
Brilliant  
56 0220

Renfert



Dynex  
Brilliant  
56 2520

Renfert



Dynex  
Brilliant  
56 0320

Renfert

## Dynex Brilliant

Flexible diamond separating discs

Item No.: 560220 - 0.2 x 20 mm (10 pk)

Item No.: 562520 - 0.25 x 20 mm (10 pk)

Item No.: 560320 - 0.3 x 20 mm (10 pk)

**BUY 4, GET 1  
FREE**

## Dynex Brilliant

Diamond grinding disc

Item No.: 560820 - 0.8 x 20 mm (1 pk)




**BUY 4, GET 1  
FREE**

# HUGE FLEXIBLE BLOCK

Flexible blocks are made from high-interlinked and high-quality nylon materials, which is practically unbreakable and flexible. Therefore, it is suitable to fabricate the framework of partial dentures, providing a more natural and comfortable options for customers.


## • Practically unbreakable and crack-resistant

HUGE Flexible block has better performance on strength compared to other nylon materials. The flexible frameworks for partial dentures are generally considered to be strong and resistant to fracture or crack.

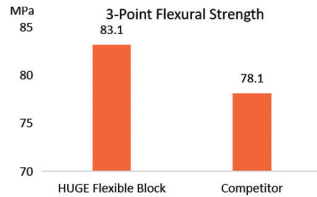


## • Stain-resistant

According to the test, nylon resin has great ability to resist stain. The sample was immersed in coffee, black tea and red wine for a period of time and no significant discoloration was observed.

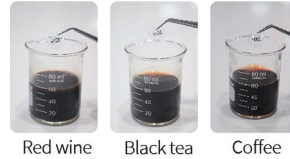


Red wine Black tea Coffee

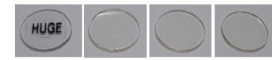

After soaking in coffee, black tea and red wine at 60 °C for 100 hours, almost no color change was detected

## • Exceptional memory to hold its form

HUGE Flexible Blocks are not easily deformable. Therefore, the framework of partial denture can maintain good elasticity and accuracy, ensuring outstanding tooth retention ability after long-term use.

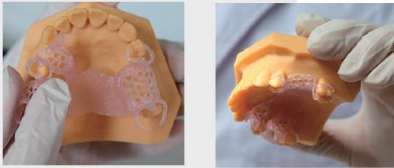

It fits good after 5,000 wearing and removal tests on model

## • Lightweight, comfortable and aesthetic

HUGE Flexible Blocks create a superior options to fabricate partial denture. The translucent and flexible materials could enhance patient wearing experience compared to metal partial dentures.




\$81.95

Available in Pink & Clear 16mm/ 20mm

# EASY TO PROCESS – THE SAME AS CONVENTIONAL PARTIAL DENTURES

## • Scan & Design

- Scan the master model using 3D scanner
- Design the framework of the partial denture from your software.

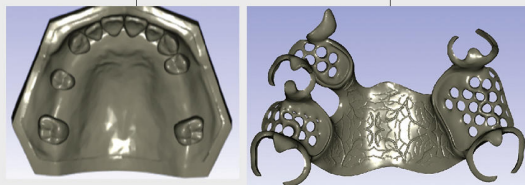

Prep Guide  
Partial Denture Framework (Minimum thickness requirements)

- Major: 1.2 mm
- Saddle Connectors Area: 0.8 mm
- Clasp on Anterior Tooth
- 3.0 mm where clasp joins framework thinning to 0.7 mm at tip
- Clasp on Posterior Tooth
- 4.0 mm where clasp joins framework thinning to 0.7 mm at tip

## • Milling


- Set the milling system according to the recommendations for PMMA.



Note:  
Preferred milling method: Wet milling  
When using dry milling method, the heat produced by the milling can cause the material to build up around the bur, Please check it regularly. If it happens, stop the milling, clean the bur and resume.

## • Add Teeth to HUGE Flexible Framework

- Arrange and wax the denture teeth on the framework.
- Process the partial denture by investing, packing and curing the acrylic resin.
- Polish the denture to ensure a smooth and comfortable fit.

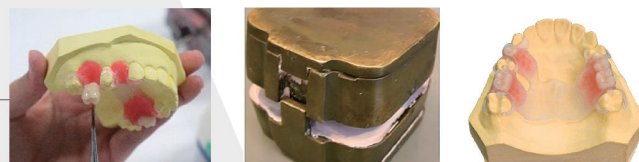

# Accufab-CEL Dental 3D Printer

**SHINING 3D**  
DENTAL

## Your Road to Success in 3D Printing

### Reliable

- Real-time cabin environment detection and heating
- The cabin temperature & humidity are continuously monitored during printing
- This heating function ensure accuracy and stability during printing, resulting in a higher success rate

### Accurate

- Accufab-CEL has achieved a high level of printing accuracy reaching  $\pm 35\mu\text{m}$ . The even distribution of light plays a crucial role in printing accuracy. 3D has developed a NEW collimated programmable lithography resulting in higher performance & accuracy.

### Affordable

- With a price of **\$7995.00** it compares favourable with the leading brands that have similar features.



## Q2 Special purchase on Accufab-CEL:

Receive 2 FREE 1 Litre bottle of Shining 3D Soft Splint Resin



**NEW**

## Introducing Soft Splint 3D Resin SS1

The optimal solution for printing soft splints with exceptional performance in Transparency, morphology and occlusion. Engineered with high tensile and flexural strength.

Great value at **\$395.00/ Litre**  
#M00005861

**SHINING 3D**  
DENTAL



# Aoralscan Elf

Smarter Occlusion

Light • Small • Powerful



Call for  
a Demo  
800.268.4442  
**\$12,600**

# TECHNICAL CORNER

## Bonding Zirconia/e-max abutments to Ti-Base using Bredent DTK-adhesive

### DTK-adhesive


- Hygienic thanks to being autoclavable
- Proper bonding thanks to FGP insulating agent
- Suitable for direct and indirect use
- Bonding of all prosthetic materials:
  - On ceramic (e.max®) → Use K-Primer
  - On zirconia/ NPM / titanium → Use MKZ Primer
  - On precious metal alloys → Use MKZ Primer + MKZ EM-Aktivator
  - On BioHPP/ PEEK → Use visio.link



### Adhesive protocol




Insulate screw channel, gingival area and shoulder of adhesive abutment and leave to dry.

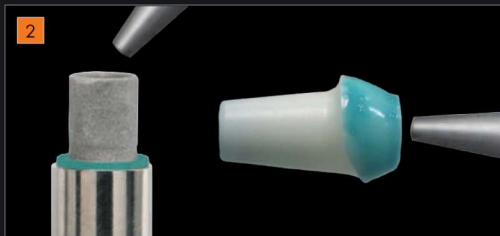

Sandblast with 110µm aluminium oxide: Pressure of 3 bar for titanium/NPM  
2 bar for zirconia



Seal screw channel with wax

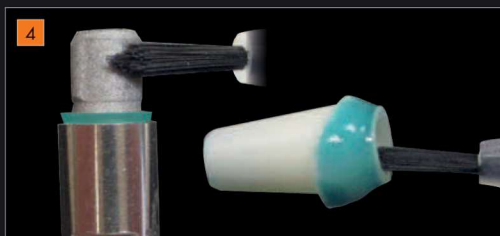


Apply MKZ Primer to titanium/NPM adhesive base and zirconia abutment.



Apply K-Primer on e-max® and MKZ Primer On titanium base.

# TECHNICAL CORNER

## Bonding Zirconia/e-max abutments to Ti-Base using Bredent DTK-adhesive

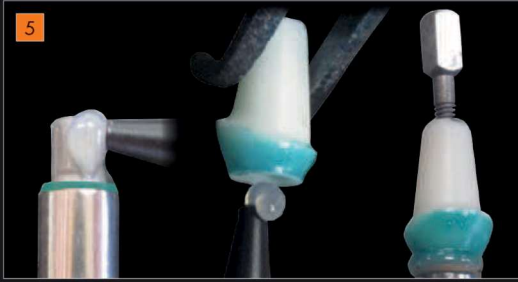

Apply DTK-adhesive (transparent or opaque), press abutment onto the adhesive base and then insert screw into the channel.



Polymerise with bre.Lux Power Unit2, for example.




Remove the screw from the channel

DTK Adhesive kit  
#54001185

Special  
Price

\$189.00  
Save  
\$30.00

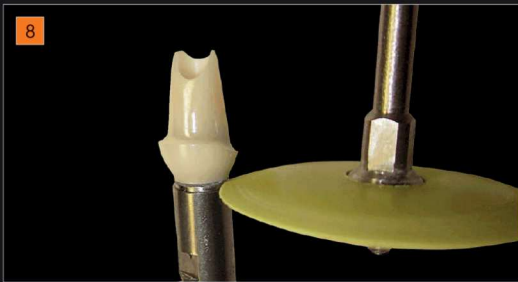

Sandblast abutment and remove excess



Bonded abutments can be sterilised.

**Note!** When using translucent and high-translucent zirconia abutments, contact manufacturer for information about hydrothermal stability achieved by aluminium oxide doping!

bredent  
group

009894GB-20180329

Mistake and subject to change reserved



All names marked with ® or ™ are protected brands and/or company names or third-party rights holders.



We take care of your smile

Ruthinium acrylic teeth  
*Premium Quality at  
an Economy Price*

Available  
exclusively

from




10+1

## AcryPlus

4 layers acrylic teeth

The nature  
in a smile.

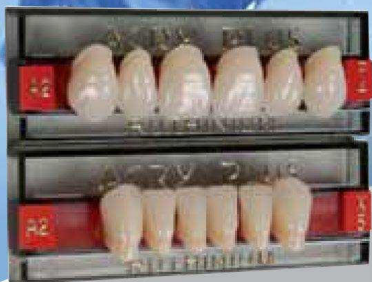

**\$28.50**

(1x6)

## AcryLux

3 layers  
acrylic teeth

Natural Bright.



**\$12.20**

(1x6)

## AcryRock

2 layers  
acrylic teeth

Natural hardness.



**\$4.99**

(1x6)

Dental Manufacturing S.p.A.

Via Cà Mignola Nuova, 1699  
45021 Badia Polesine - Rovigo (Italy)  
+39 0425 51628 - 0425 590587 - Fax +39 0425 590156  
info@ruthinium.it

[www.ruthinium.it](http://www.ruthinium.it)



# Flexing to meet every shade needs

## PMMA Discs (Shaded and Multi-layer)



**2+1**

BUY 2 PMMA Discs,  
**GET 1 FREE.**  
Promo code: CPMMAFY26

Dentsply Sirona offer practices and laboratories a continuously growing variety of materials for economical and precise fabrication of clinically flawless and esthetically high-quality prosthetics.

## Cercon® Discs



**3+1**

BUY 3 Cercon Discs,  
**GET 1 FREE.**  
Promo code: CCER31FY26

**20+5**

BUY 20 Cercon Discs,  
**GET 5 FREE.**  
Promo code: CCER20FY26

## Lucitone Digital Fit™ Discs



**3+1**

BUY 3 Lucitone Digital Fit Discs,  
**GET 1 FREE.**  
Promo code: CLFD31FY26

**20+5**

BUY 20 Lucitone Digital Fit Discs,  
**GET 5 FREE.**  
Promo code: CLFD20FY26



**PRECISE PASTE**  
 WR600 - \$53.99/ 3oz paste  
 WR602 - \$169.99/ 12oz paste  
 PRECISE Pressure Indicator Paste is used for seating full, immediate and partial dentures.

**PRECISE BRUSHES**  
 WR604 - \$17.50  
 Aluminum Applicator Brushes – The Perfect Partner for PRECISE Pressure Indicator Paste.



**PRECISE Brush Cleaner**  
 WR603 - \$23.45/ 150 ml Bottle  
 Quickly dissolves any type of pressure indicator paste and wipes away clean.



**SMOOTH-EX**  
 WR750 - \$23.25  
 The excellent solution for creating the perfect finished edge on mouthguard and bleaching materials.



**Round-Off**  
 WR760 - \$22.95  
 The essential solution for rounding off edges of dual laminate materials and hard vacuum forming resin.



**E-Z-D Flasking Sheets**  
 WR980 - \$14.98  
 Pre-cut paper flasking sheets are valuable time savers. Sheets allow for easy deflasking and allow easy identification of the cases before complete deflasking is done.



**Stik-Bak**  
 #0705 - \$15.49  
 Specially designed for remounting denture models - Stik-Bak offers an incredibly strong hold with just a few drops and will separate from the mounting with a light tap from a plaster knife and hammer.



Isonic Ultrasonic Cleaner  
#DS400B

Includes 600 ml Beaker

**SPECIAL PRICE - \$195.00**



Sunburst MicroMotor  
#P-900/ X-99L

- Light Weight
- Great for Chairside Adjustments

**SPECIAL PRICE - \$350.50**



Vibrator  
#307

**SPECIAL PRICE - \$125.00**



Meta Grinding Box  
#800-2100

**SPECIAL PRICE - \$349.00**



Angle Plier  
#070351 -  
**SPECIAL PRICE - \$65.50**

3 prong Plier  
#070352 -  
**SPECIAL PRICE - \$65.00**



Plaster Nippers  
#650-19

**SPECIAL PRICE - \$85.00**



Acrylic Micing Bowl Large with Lid  
#D11A

**SPECIAL PRICE - \$35.00**



Mestra Denture Duplicating Flask  
#020206

**SPECIAL PRICE - \$195.00**

**2+1**



Alginate Mixer  
#ADA004002

**SPECIAL PRICE - \$1495.00**

**SAVE \$300**



Mestra Pour Flask  
#020202

**SPECIAL PRICE - \$175.00**



Nitrile - 3.2 mil  
Ice Blue

**8+2**

Nitrile - 2.5 mil  
Aqua

**8+2**

Latex - 7.2 g  
Creamy White

**8+2**

- XS 9889A5
- S 9889A6
- M 9889A7
- L 9889A8
- XL 9889A9

**100 per Box**

- XS 92395
- S 92396
- M 92397
- L 92398
- XL 92399

**200 per Box**

- XS 29225
- S 29226
- M 29227
- L 29228
- XL 29229

**100 per Box**

Price: \$10.<sup>00</sup>  
Net Price: \$8.<sup>00</sup>

Price: \$19.<sup>65</sup>  
Net Price: \$15.<sup>70</sup>

Price: \$13.<sup>00</sup>  
Net Price: \$10.<sup>40</sup>





VITA



**VITA VIONIC® BOND multiDose COMPLETE**  
Limited Time Offer


Strengthen your digital denture workflow with a bonding solution you can truly rely on. **VITA VIONIC BOND multiDose COMPLETE** delivers **consistent, high-performance results** designed to support precision, longevity, and efficiency in every case. **Elevate your workflow, boost your productivity.** Take advantage of this exceptional offer by contacting Central Dental today



Measure


Apply

Bond

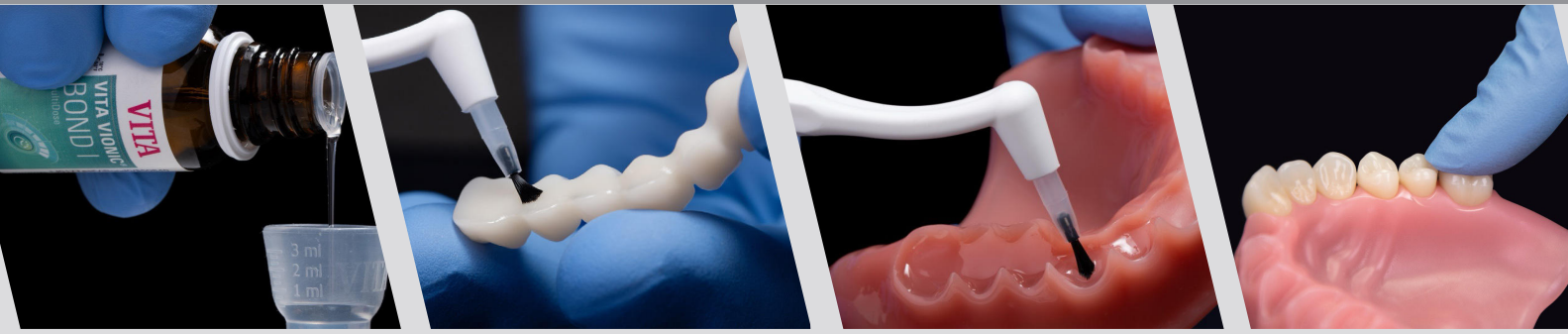

Buy 10 **VITA VIONIC® DENT DISCS**,  
Get a free **VITA VIONIC BOND multiDose!**

Promo Code: QDEN2622

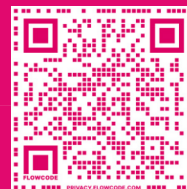

Scan the QR code to  
watch a short tutorial  
on a **VITA VIONIC**  
Workflow!

© 2026 VITA North America. Valid on purchases made April 1-June 30, 2026 only. Offers valid in Canada only. Offers may not be combined with any other promotion or discounts. Invoice must be submitted to [promo@vitanorthamerica.com](mailto:promo@vitanorthamerica.com) no later than July 15, 2026. VITA North America reserves the right to discontinue or change the offers at any time without notice. No cash value.

**800-828-3839**  
[www.vitanorthamerica.com](http://www.vitanorthamerica.com)



# NU-DENT<sup>®</sup>

## PROFESSIONAL STRENGTH DENTURE & APPLIANCE CLEANER

Get 4  
Bottles  
FREE

with  
purchase  
of a case  
(48)

For the safe and effective  
removal of coffee, tea, tobacco,  
tartar and calculus  
without damaging  
the appliance.



**64 Cleanings per bottle**

Cleans in Less than 25 min.



**CD** *Central Dental*  
The laboratory and dentist supply specialist.

Tel: 1-800-268-4442 Fax: 416.694.1071 [www.CentralDentalLtd.com](http://www.CentralDentalLtd.com)

Not recommended for use on Gold.